

High-Performance Thin-Layer Chromatography: A Fast and Efficient Fingerprint Analysis Method for Medicinal Plants

HPTLC Fingerprint of Ginsenosides

Tiên Do, Scientific Support Specialist / Deputy head of laboratory (CAMAG), tien.do@camag.com

Melanie Broszat, Scientific Business Development Manager (CAMAG), melanie.broszat@camag.com

Matthias Nold, Product Manager Reference Materials, matthias.nold@sial.com



Ginsenosides are triterpene saponins. Most ginsenosides are composed of a dammarane skeleton (17 carbons in a four-ring structure) with various sugar moieties (e.g. glucose, rhamnose, xylose and arabinose) attached to the C-3 and C-20 positions.

Over 30 ginsenosides have been identified and classified into two categories:

- 20(S)- protopanaxadiol (PPD) (Rb1, Rb2, Rb3, Rc, Rd, Rg3, Rh2, Rs1)
- 20(S)-protopanaxatriol (PPT) (Re, Rf, Rg1, Rg2, Rh1). The difference between PPTs and PPDs is the presence of a carboxyl group at the C-6 position in PPDs.

Moreover, several rare ginsenosides, such as the ocotillol saponin F11 (24-R-pseudoginsenoside) and the pentacyclic oleanane saponin Ro (3,28-o-bisdesmoside) have also been identified. Asian ginseng (*Panax ginseng*) is commonly known as the true ginseng.

In this article, HPTLC methods suitable for the analysis of ginsenosides are presented, using CAMAG equipment, Merck TLC plates, analytical standards and extract reference materials. The extract reference materials are manufactured by HWI Analytik and exclusively distributed by Merck Sigma-Aldrich.

Detection of ginsenosides in the HPTLC fingerprint of different *Panax* species (roots and root extracts) is obtained by following the HPTLC methods of the Ph. Eur. monograph,¹ USP DSC 2015 monograph² and

the method of the HPTLC Association (International Association for the Advancement of HPTLC)³ by comparison of the R_f values and colors of reference substances and matching zones in the root extract. Depending on the regulation followed, one of the three methods of identification can be used.

Recommended CAMAG Devices:

Automatic TLC Sampler (ATS 4), Automatic Developing Chamber (ADC 2), TLC Visualizer 2, Chromatogram Immersion Device 3, TLC Plate Heater 3, and visionCATS

Derivatization Reagent:

Anisaldehyde¹ or sulfuric acid reagent^{2,3}

Sample:

0.015 g/mL extract (HWI extract) in 70% methanol¹

Note: Deviation from methods 2 and 3 for the sample preparation

Standards:

Standard solutions of ginsenosides were prepared in a concentration of 0.5 mg/mL in methanol.

Note: Deviation from method 3 for the application volumes

Chromatography Following USP <203> ⁴:

Stationary phase: HPTLC Si 60 F₂₅₄ 20 x 10 cm (Merck)

Sample application: 4 μ L each of test solution and 2 μ L of standards are applied as 8 mm bands, 8 mm from lower edge, 20 mm from the left edge^{1, 3}

Note: Deviation from method 3 for the sample preparation and the application volume of reference and test solutions

4 μ L each of test solution and 4 μ L of standards are applied as 8 mm bands, 8 mm from lower edge, 20 mm from the left edge²

Note: Deviation from method 2 for the sample preparation

Developing solvent: Ethyl acetate – water – butanol
25:50:100 (v/v/v) - upper layer¹

Dichloromethane – ethanol – water 70:45:6.5 (v/v/v)²

Chloroform – ethyl acetate – methanol – water
15:40:22:9 (v/v/v/v)³

Development: In the ADC 2 with an unsaturated chamber and after conditioning at 33% relative humidity for 10 min using a saturated solution of magnesium chloride¹

Development is performed with ADC 2, saturated for 20 minutes with the developing solvent (filter paper). Prior to the development the plate is conditioned for 10 min to a relative humidity of 33% (with a saturated solution of $MgCl_2$).^{2, 3}

Developing distance: 70 mm (from lower edge)^{1, 2}
80 mm (from lower edge)³

Plate drying: 5 min in a stream of cold air

Derivatization: The plate is immersed (immersion speed: 3 cm/s, immersion time: 0 s) into anisaldehyde reagent (mixture of 0.5 mL of *p*-anisaldehyde, 10 mL of glacial acetic acid, 85 mL of methanol, and 5 mL of sulfuric acid) with the Chromatogram Immersion Device 3 and heated for 5 min at 105°C.¹

The plate is immersed (immersion speed: 3 cm/s, immersion time: 0 s) into sulfuric acid reagent (10% in methanol) with the Chromatogram Immersion Device 3 and heated for 5 min at 100 °C.^{2, 3}


Evaluation: Documentation under method 3 is with 100 °C white light^{1, 2, 3} and UV 366 nm^{2, 3} after derivatization with the TLC Visualizer 2

Results and Discussion:

HPTLC chromatograms of ginsenoside standards and a *Panax ginseng* root extract


HPTLC chromatograms after derivatization. Tracks 1-17: ginsenosides, track 18: protopanaxadiol, track 19: *Panax ginseng* root extract (article no.: 05115001 batch: HWI01294)

Method 1: According to the Ph. Eur.¹

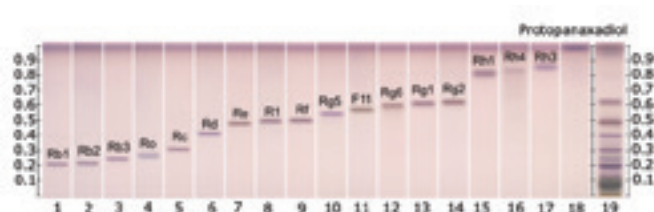


HPTLC chromatograms under white light after derivatization

Method 2: According to the Tianshi ginseng method from²

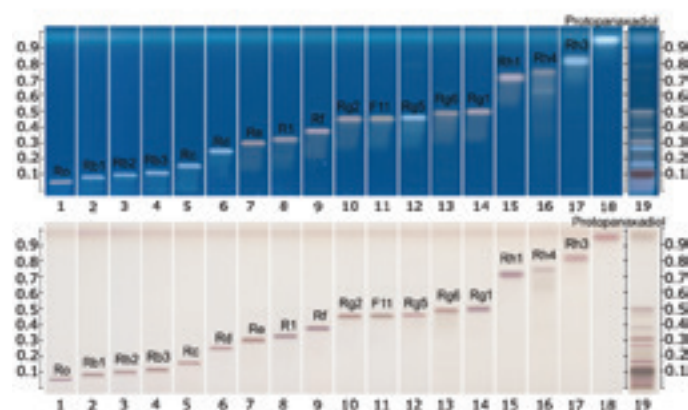

HPTLC chromatograms under UV 366 nm and under white light after derivatization

Method 3: According to the HPTLC Association³

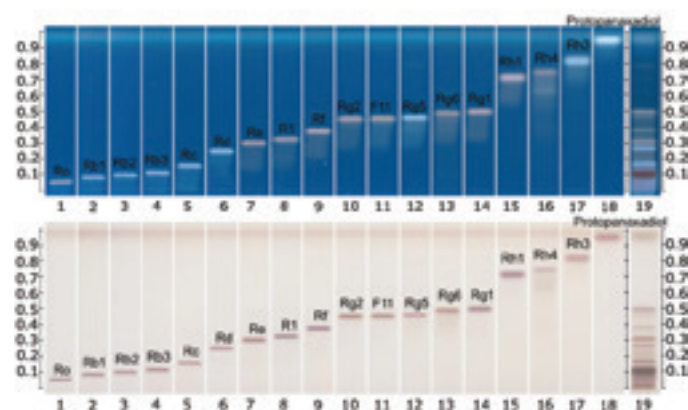

HPTLC chromatograms under UV 366 nm and under white light after derivatization.

HPTLC chromatograms of plants containing ginsenosides (different *Panax* species)

P. ginseng, *P. quinquefolium*, *P. notoginseng*, *P. japonicus*, *P. vietnamensis* roots and root extracts were collected and analyzed. The *P. ginseng* root extract (article no.: 0511-50-01 batch: HWI01294) was used as botanical reference material to identify Asian ginseng (the ginsenoside Rf should be present and F11 absent).


Did you know...

...that Merck also offers dedicated MS-grade TLC & HPTLC plates for TLC-MS coupling?

Contact us to request your copy of our Brochure "TLC and HPTLC MS-grade plates"

[SigmaAldrich.com/tlc](https://sigmaaldrich.com/tlc)


Method 1: According to the Ph. Eur.¹

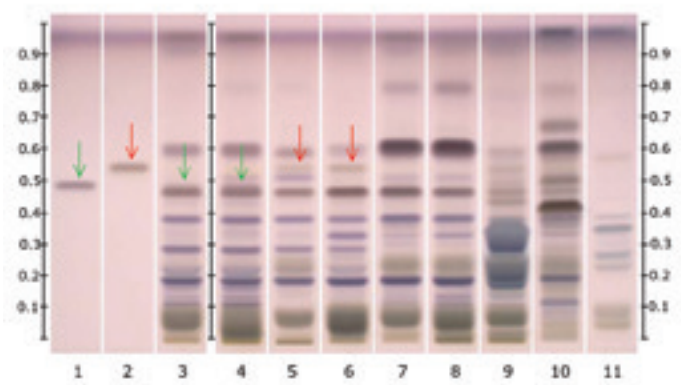


HPTLC chromatogram under white light after derivatization.

Track 1: ginsenoside Rf (green arrow); **2:** ginsenoside F11 (red arrow); **3:** *P. ginseng* root extract (article no.: 0511-50-01 batch: HWI01294, ginsenoside Rf highlighted with the green arrow); **4:** *P. ginseng* root (ginsenoside Rf highlighted with the green arrow); **5:** *P. quinquefolium* root extract (American ginseng, ginsenoside F11 highlighted with the red arrow); **6:** *P. quinquefolium* root (American ginseng, ginsenoside F11 highlighted with the red arrow); **7:** *P. notoginseng* root extract; **8:** *P. notoginseng* root; **9:** *P. japonicas* root; **10:** *P. vietnamensis* root; **11:** wild Vietnamese ginseng root.

Method 2: According to the Tianshi ginseng method from²

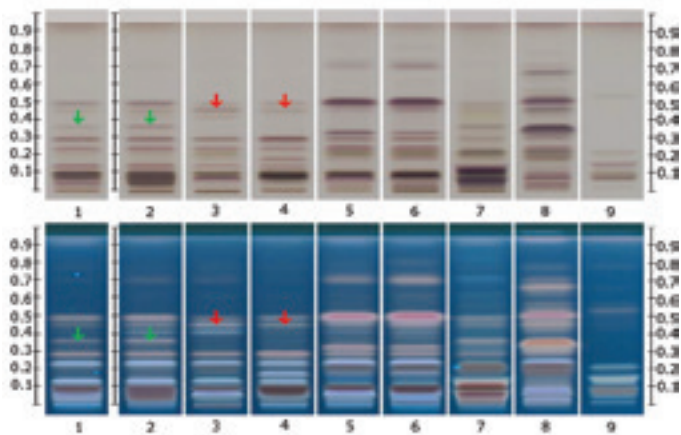

Method 3: According to the HPTLC Association³

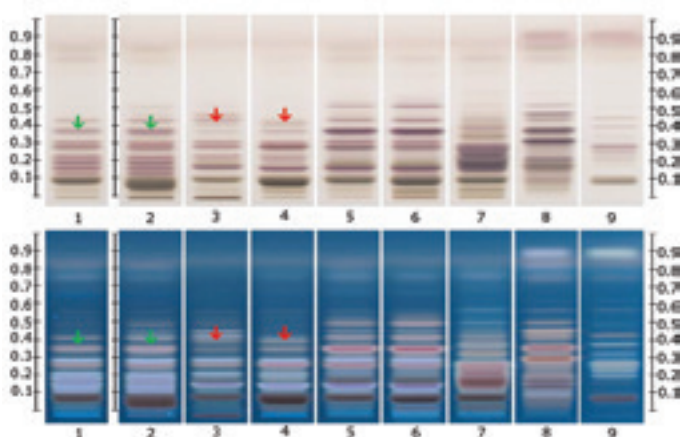

Method 2 & 3 Parameter:

HPTLC chromatograms under white light and UV 366 nm after derivatization with ginsenoside Rf (green arrows) and ginsenoside F11 (red arrows).

Track 1: *P. ginseng* root extract (article no.: 05115001 batch: HWI01294, ginsenoside Rf highlighted with the green arrows); **2:** *P. ginseng* root (ginsenoside Rf highlighted with the green arrows); **3:** *P. quinquefolium* root extract (American ginseng, ginsenoside F11 highlighted with the red arrows); **4:** *P. quinquefolium* root (American ginseng, ginsenoside F11 highlighted with the red arrows); **5:** *P. notoginseng* root extract; **6:** *P. notoginseng* root; **7:** *P. japonicas* root; **8:** *P. vietnamensis* root; **9:** wild Vietnamese ginseng root.

All shown methodologies are suitable for detection of ginsenosides in different *Panax* species. Ginsenoside Rf is unique to Asian ginseng while F11 is found exclusively in American ginseng. Thus the Rf/F11 ratio is used as a phytochemical marker to distinguish American ginseng from Asian ginseng. In the botanical reference material used (*Panax ginseng* root extract, article no.: 0511-50-01) the presence of Rf and absence of F11 could be confirmed with all methods and it is therefore suited for identification of Asian ginseng.

Methods 1 and 2 have the advantage in that they are in accordance with official monographs in the pharmacopoeias (Ph. Eur. resp. USP). Method 3 is an alternative method provided by the HPTLC Association. The method is improved for the separation of Rf and F11 to better distinguish American and Asian ginseng. The derivatization with sulfuric acid reagent (methods 2 and 3) leads to different colored zones under UV 366 nm, useful for identification.

Featured Products

Description	Package Size	Cat. No.
Analytical Standards		
Ginsenoside Rb1*	10mg	170580
Ginsenoside Rb2	10mg	41868
Ginsenoside Rb3	10mg	42635
Ginsenoside Rc	5mg	44987
Ginsenoside Rd	10mg	01518
Ginsenoside Re*	10mg	03000590
Ginsenoside Re	10mg	77960
Ginsenoside Rf*	10mg	01580590
Ginsenoside Rg1*	10mg	00370580
Ginsenoside Rg2	10mg	08171
Ginsenoside Rg3	10mg	64139
Ginsenoside Rg5	5mg	43016
Ginsenoside Rh1	10mg	56805
Ginsenoside Rh3	5mg	43084
Ginsenoside Rh4	5mg	42776
Ginsenoside Ro	10mg	94381
Notoginsenoside R1	10mg	77089
Protopanaxadiol	10mg	62685
Protopanaxatriol	10mg	42476
Pseudoginsenoside F11	10mg	67530

*HWI reference standard

(continued on next page)

Featured Products (cont.)

Description	Package Size	Cat. No.
Extract Reference Material		
Panax ginseng extract	500mg	05115001
Quantitative Markers: Ginsenoside Rb1, Rb2, Rc, Rd, Re, Rf, Rg1 and Rg2		
Qualitative Markers**: Ginsenoside Rb1 and Rg1		
(**: traceable to HWI primary pharmaceutical standards)		
TLC Plates		
HPTLC glass plate Silica gel 60 F254 20×10 cm	50 Plates	1.05642

References

1. Ginseng root: Monograph in Ph. Eur. 8.0, 01/2008:2383. European Directorate for the Quality of Medicines and HealthCare, Strasbourg, France.
2. Tienchi Ginseng Root and Rhizome Dry Extract: Monograph in USP 40-NF35 (2017). United States Pharmacopeial Convention, Rockville, MD, USA.
3. HPTLC identification method for *Panax ginseng*, HPTLC Association (www.hptlc-association.org) (accessed May 22, 2016)
4. <203> High-Performance Thin-Layer Chromatography Procedure for Identification of Articles of Botanical Origin in USP 39-NF34. (2016), United States Pharmacopeial Convention, Rockville, MD, USA

Find all available analytical standards for phytochemicals listed on

SigmaAldrich.com/medicinalplants

For all ginsenoside standards see

SigmaAldrich.com/panax

An overview of all plant extract reference materials can be found at SigmaAldrich.com/plantextracts

NUTRITIONAL SUPPLEMENTS

New Phytochemical Standards from HWI Pharma Solutions

Matthias Nold, Product Manager, Reference Materials, matthias.nold@sial.com



We are proud to be an exclusive distributor of the reference standards manufactured by HWI-Pharma Solutions in Rülzheim (Germany).¹ Since the launch of the first series of products in January 2011, the product range has been continuously expanded and currently consists of more than 120 products.

The quantitative value is determined by **quantitative NMR** (qNMR).² This is a direct relative method, which is increasingly used for the quantification of organic compounds, as an alternative to the much more laborious mass balance approach. The certificate

delivered with these products also contains a chromatographic purity value.

Recently, several new products have been added to the portfolio (see **Table 1**). A complete listing can be found online at SigmaAldrich.com/phytopharma.

Table 1. Newly added HWI reference standards of phytopharmaceuticals

Description	Package Size	Cat. No.
Benzyl acetate	100mg	05880595
Berberine chloride	50mg	00900585
Carminic acid	25mg	03320585
(±)-β-Citronellol	100mg	05630590
Ectoïne	100mg	02380595
Hydroxyectoïne	100mg	02390595
Isoxanthohumol	25mg	05890580
DL-Kavain	25mg	05790585
Alpha-onocerin	10mg	05800590
Patchouli alcohol	10mg	05690595

References:

1. G. Förster, F. Michel, M. Nold; Analytix 1/2010 page 11.
2. M. Veith; Analytix 1/2010 page 14.